

# Magnetic Ear Tag Reporter for Automated Quantification of Head Twitch Response – Setup Manual

Kwan Lab, Yale University



<http://alexkwanlab.org/>

Last updated: July 27, 2021

## Table of Contents

<b>Overview .....</b>	<b>2</b>
<b>Section 1: Box, Coil, and Phono Preamp.....</b>	<b>4</b>
Box, Coil, Phono Preamp, and Wiring .....	4
Parts List.....	6
Assembly .....	6
<b>Section 2: Data Acquisition Device (DAQ).....</b>	<b>6</b>
DAQ and Wiring.....	6
Parts List.....	7
Assembly .....	7
Notes .....	8
<b>Section 3: Ear-tagging the Mouse .....</b>	<b>8</b>
Magnetic Ear-tagging.....	8
Parts List.....	8
Assembly .....	9
Notes .....	9
<b>Section 4: Optional Video Recording Setup.....</b>	<b>9</b>
Recording Setup.....	9
Parts List.....	10
Assembly .....	11
Notes .....	11


<b>Section 4: Computer Control.....</b>	<b>11</b>
Control of the DAQ and Camera via MATLAB .....	11
Parts List.....	11
Assembly .....	11
Notes .....	12
<b>Addendum: Scaling Up .....</b>	<b>12</b>
Changes to the Setup when Scaling Up .....	12

## Overview


This document provides instructions for creating a setup to automatically quantify the head twitch response in mice via a magnetic ear tag. When administered psilocybin, LSD, DOI, or most other hallucinogens, mice exhibit a characteristic head twitch caused primarily by the binding of the drug to the serotonin 5-HT<sub>2A</sub> receptor. Previously, the only way to quantify these head twitches was to videotape the mouse's behavior, which could then be visually analyzed by a trained observer; this method is tedious and time-consuming.

A new, simpler, and minimally invasive method of quantifying the head twitch response was described by Revenga et al. in 2020 and relies on the principle of electromagnetic induction. From Faraday's Law, a changing magnetic flux (which is a measure of the magnetic field passing through an area) in a closed connecting loop of coil induces a momentary current within the coil. One way to change the magnetic flux in order to generate the induced current is by quickly shaking a magnet within the loop of coil. The induced electromotive force—i.e., the voltage source—responsible for the induced current to flow can then be quantified using a voltmeter or a data acquisition system. When a magnetic ear tag is attached to a mouse's ear and the mouse twitches its head upon administration of the drug, the rapid shaking of the magnet leads to a characteristic voltage spike that can be identified via a MATLAB script. This document will describe how to build a magnetic ear tag reporter setup, based on the 2020 Revenga et al. paper, to automatically quantify this head twitch response.

A completed setup looks like this:

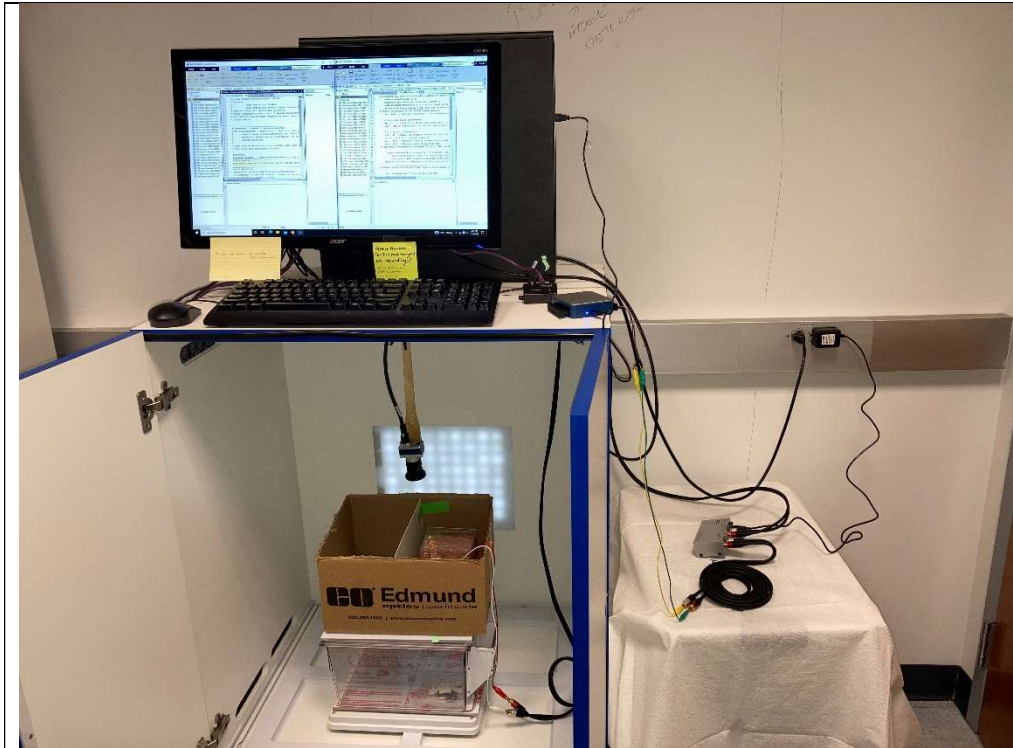


*A view of the entire setup*

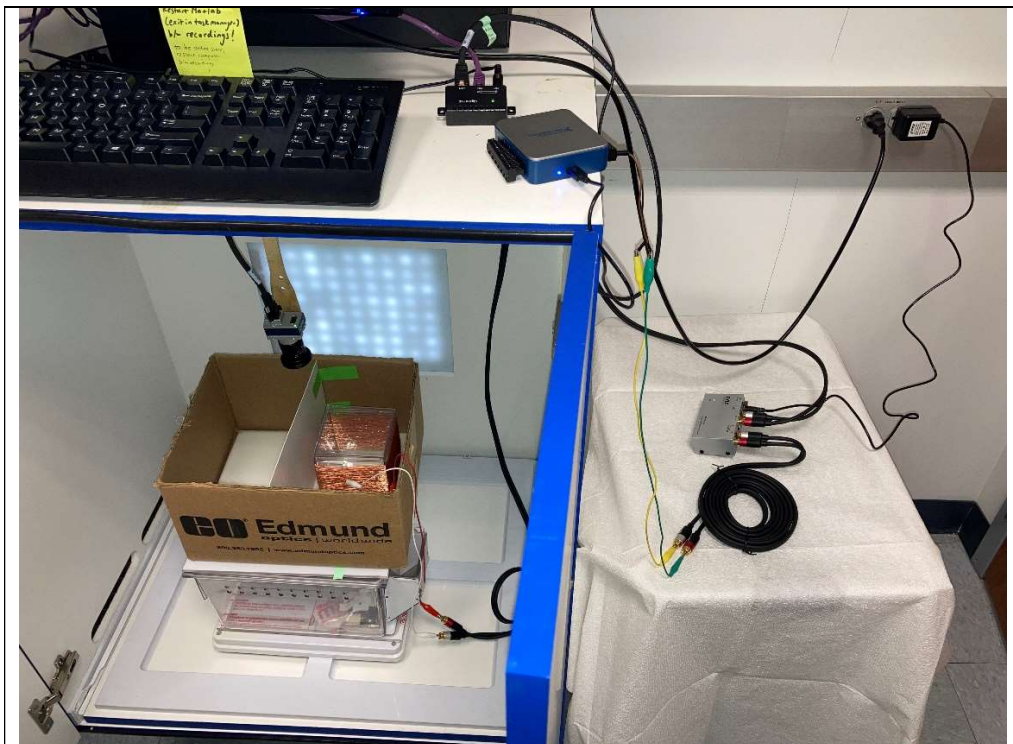

*A closeup view of the box & coil attached to the phono preamp and DAQ. The wires connecting the box to the phono preamp and the wire connecting the camera to the power injector go through a hole in the top right of the crate.*


Each section in this document describes a component of the system, and the diagram below illustrates how these components are connected together:

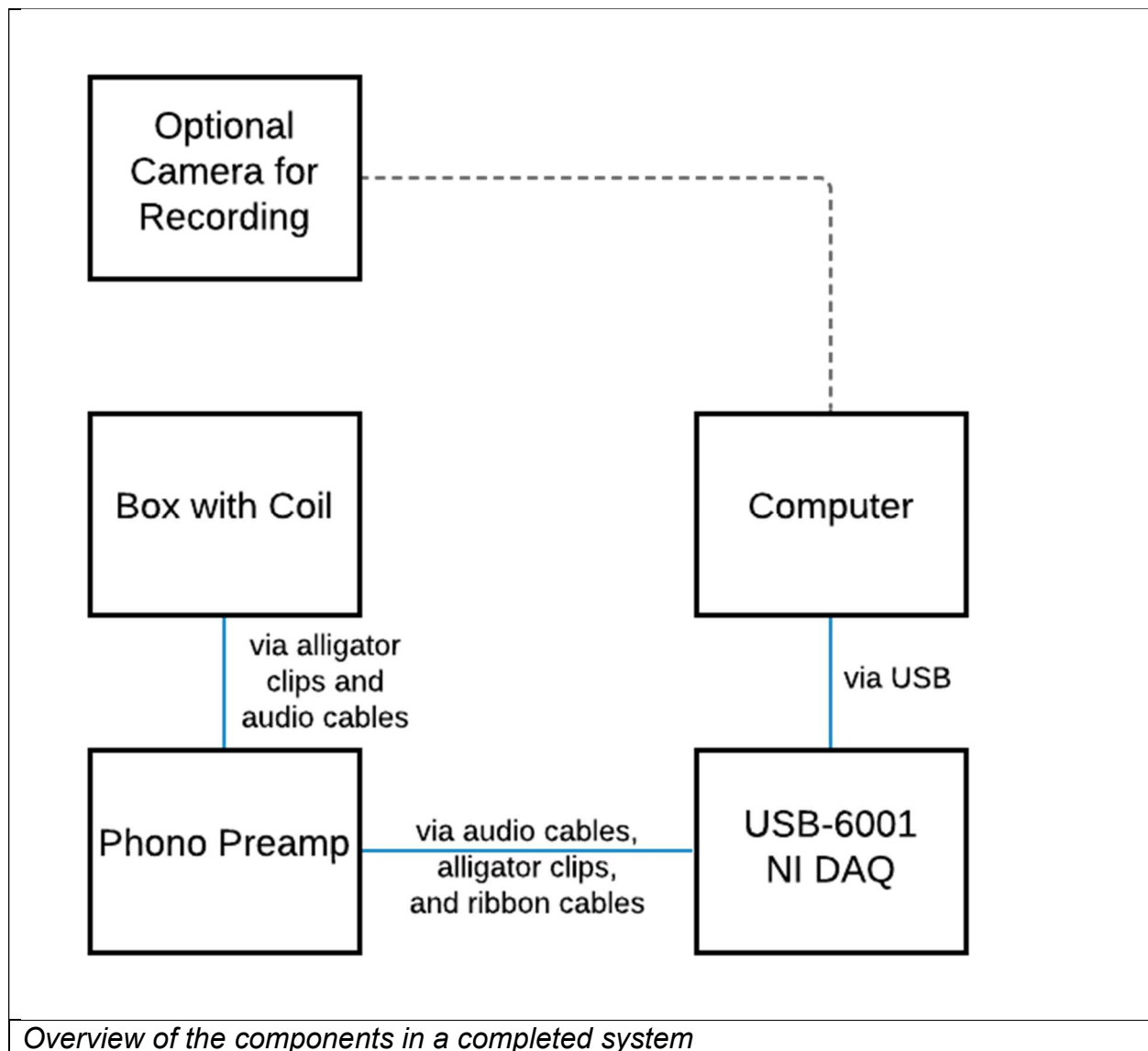


## Section 1: Box, Coil, and Phono Preamp

### Box, Coil, Phono Preamp, and Wiring

Following intraperitoneal injection with the drug, the mouse with the magnetic ear tag will roam around in the plastic box wrapped around with coil. The phono preamp amplifies the voltage signal.

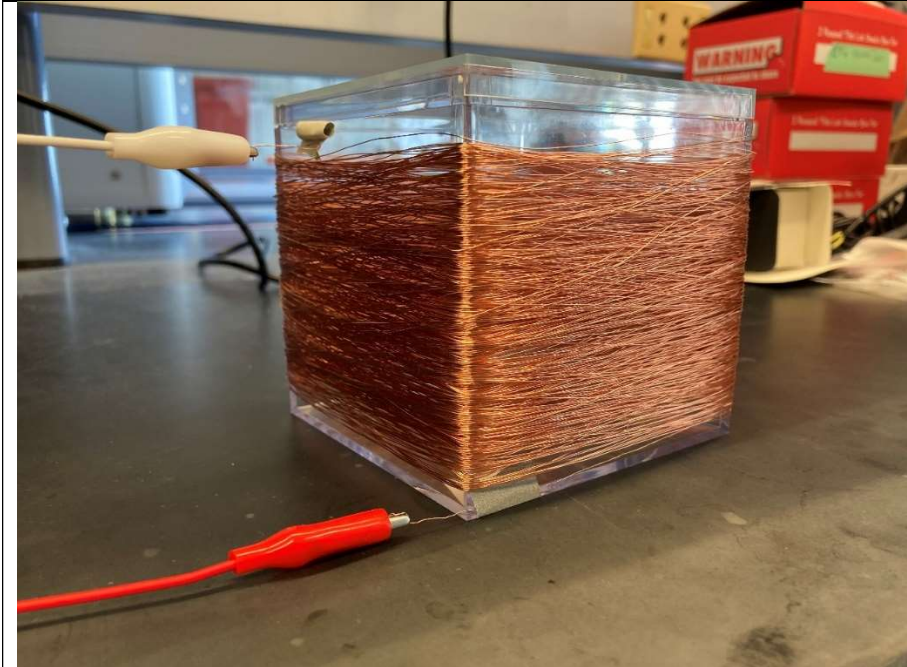

*Closeup of the plastic box wrapped with coil, the ends of which are attached to alligator clips*

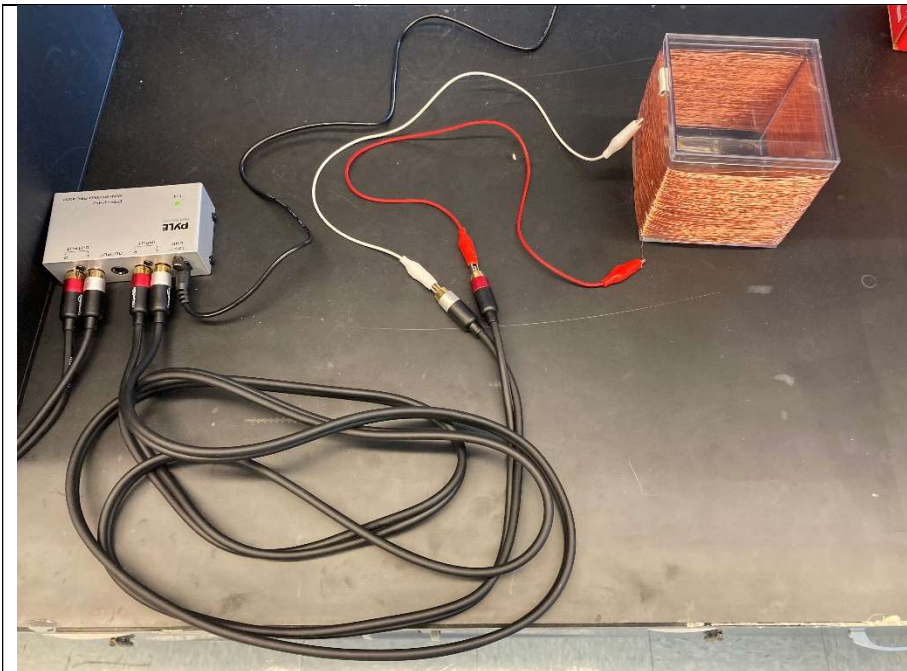

*Connecting the wires around the plastic box to the phono amp input using audio cables*



## Parts List

- Clear plastic box cube (4" x 4" x 4")
  - <https://www.amazon.com/Clear-Plastic-Box-Cube-Square/dp/B0883FNSFJ/>
- 30 AWG copper wire
  - <https://www.amazon.com/BNTECHGO-AWG-Magnet-Wire-Transformers/dp/B00UWCXRK6>
- Phono preamp
  - <https://www.amazon.com/Pyle-Phono-Turntable-Preamp-Preamplifier/dp/B004HJ1TTQ>
- Two wires with alligator clip on each end
  - <https://www.amazon.com/WGGE-WG-026-Pieces-Colors-Alligator/dp/B06XX25HFX>
- Set of 2-male to 2-male audio cables
  - <https://www.amazon.com/AmazonBasics-2-Male-RCA-Audio-Cable/dp/B01D5H8P0G>


## Assembly

1. Wrap all ~820 feet of the copper wire around the clear plastic box, ensuring that both ends of the wire are easily accessible.
2. Connect the alligator clips to the two ends of the wire. Attach the alligator clips on the other side of the wires to the ends of the audio cables.
3. Plug the audio cables into the phono preamp, which should be plugged into the wall socket.

## Section 2: Data Acquisition Device (DAQ)

### DAQ and Wiring

A data acquisition device (DAQ) is used to track the voltage signal over time.

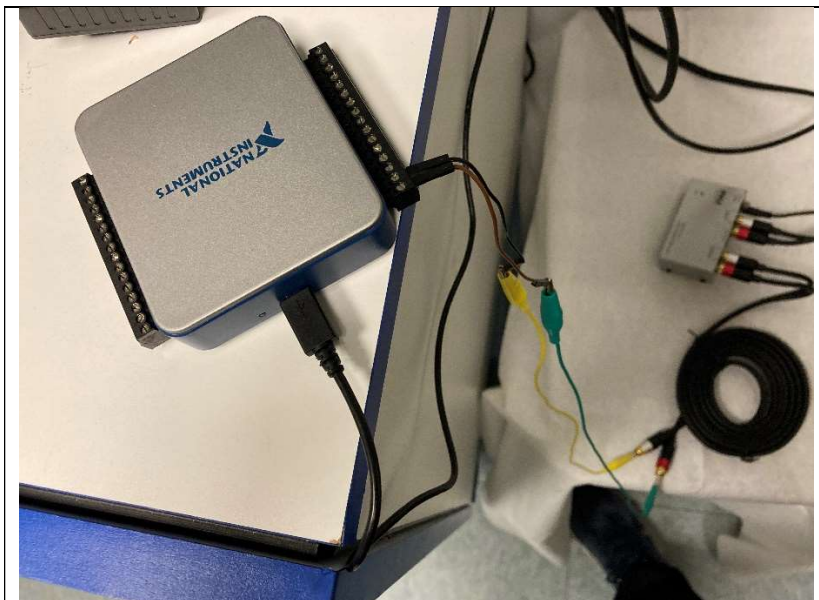

*DAQ attached via ribbon cables and alligator cables to the phono preamp. A USB cable connects to the computer.*



*Closeup of the DAQ with ribbon cables*

#### Parts List

- NI USB-6001 DAQ
  - <https://www.ni.com/en-us/shop/hardware/products/multifunction-io-device.html?modelId=124894>
- Flathead screwdriver
- Two wires with alligator clip on each end
  - <https://www.amazon.com/WGGE-WG-026-Pieces-Colors-Alligator/dp/B06XX25HFX>
- Set of 2-male to 2-male audio cables
  - <https://www.amazon.com/AmazonBasics-2-Male-RCA-Audio-Cable/dp/B01D5H8P0G>
- Male to male ribbon cables
  - <https://www.amazon.com/Elegoo-EL-CP-004-Multicolored-Breadboard-arduino/dp/B01EV70C78/>

#### Assembly

1. Plug the USB cable from the DAQ into the computer.
2. Connect the audio cables and alligator cables to the ribbon cables.
3. Use a flathead screwdriver to insert the ends of the ribbon cable into the analog input (AI) channels. In the image above, the ribbon cables are attached to the analog input channel 0 (in the + and – ports).
4. Use the MATLAB script described in Section 4 to allow the DAQ to record voltage signals.

## Notes

- If the voltage signal from the DAQ appears to stay the same and does not have any peaks, make sure that the phono amp is plugged in and attached to the DAQ. The voltage signal from a head twitch is not recognizable without the phono preamp.
- While other DAQ's may work, the NI DAQ is especially helpful because it can be directly controlled via MATLAB.

## Section 3: Ear-tagging the Mouse

### Magnetic Ear-tagging

Each mouse will have an ear tag superglued to a magnet attached to one of their ears.



*Ear tag with magnet attached to a mouse's ear*

## Parts List

- La Pias Aluminum Ear Tags (Item 56780-100 from Stoelting)
  - <https://stoeltingco.com/Neuroscience/Ear-Tags--Mouse~9785>
- Ear Tag Applicator (Item 56791 from Stoelting)
  - <https://stoeltingco.com/Neuroscience/Ear-Tags--Mouse~9785>
- Neodymium Magnets (Product D1005-10 from Super Magnet Man; N45 magnet with diameter 3mm and thickness 0.5mm)
  - [https://supermagnetman.com/products/d1005-10?\\_pos=1&\\_sid=55a9637bc&\\_ss=r](https://supermagnetman.com/products/d1005-10?_pos=1&_sid=55a9637bc&_ss=r)
- Superglue
  - <https://www.amazon.com/Loctite-Control-4-Gram-Bottle-1739050/dp/B00ELV2D0Y>
- Nitrocellulose Marker



- <https://www.amazon.com/Colortone-Touch-up-Marker-Opaque-Lacquer/dp/B0731FS9YS>
- Plastic Forceps
  - <https://www.finescience.com/en-US/Products/Forceps-Hemostats/Standard-Forceps/Plastic-Forceps/11700-00>

### Assembly

1. Use the superglue to attach the male end of the ear tag to the magnet. Coat the magnet with the nitrocellulose marker (in order to prevent irritation to the mouse's ear caused by the magnet). Wait a few hours to dry.
2. When it is time to ear tag the mouse, place the female end of the ear tag onto the red slot of the applicator and the male end of the ear tag with the magnet on the other slot. Briefly place the mouse under anesthesia (via isoflurane) and apply force to the handles of the applicator to position the ear tag in the middle part of the mouse's ear (avoiding the cartilage) such that the magnet faces inward.
3. Allow around 2 days for the animal to rest before experimenting (to prevent excessive ear grooming).

### Notes


- The neodymium N45 magnet (product D1005-10 from Super Magnet Man with a diameter of 3mm and a thickness of 0.5mm) was carefully selected based on trial and error. Initial trials used the neodymium N50 magnet (product D1007 from Super Magnet Man with a diameter of 3 mm and a thickness of 1 mm), but when magnets were attached to both ears of one mouse, the two ears stuck together. Furthermore, when the magnet was attached to one ear of each mouse within group housing, the mice would often get magnetized together, preventing them from moving independently and irritating their ears; sometimes the mice would pull each other's ear tags off completely. The voltage signal from the thinner, N45 magnet was just as prominent as the thicker magnet and the mice rarely got magnetized to one another; when they did, they were easily able to separate.
- When mice are group-housed, ensure that the same pole on each magnet faces the outward direction to prevent the mice from getting magnetized together.
- It is not necessary to attach an ear tag with the magnet to both of a mouse's ears. One ear tag will suffice, and it not a good idea to have two ear tags as that could cause the mouse's ears to be magnetized to one another.
- After supergluing the magnet to the ear tag, wait a few hours before applying the ear tag on the mice. If the glue is not dry enough, the magnet might get stuck to the magnetic applicator instead of the ear tag. Furthermore, after applying the ear tag to the mouse, use plastic forceps (instead of magnetic ones) to separate the magnet on the ear tag from the applicator.
- If the mouse is not given enough time to rest before experimenting, it might excessively groom its ear, leading to potential false positives.

## Section 4: Optional Video Recording Setup

### Recording Setup

Ideally, a camera will not be necessary because the magnetic ear tag setup can quantify head twitches. However, to verify that the setup works or in order to add a second layer of validity to

the quantification of the head twitch response, it may be useful to incorporate a video recording device. A camera able to record a high frame rate (around 180 fps) is necessary to visualize the high-frequency head twitches. We use a Dalsa Genie Nano GigE PoE camera connected to a power injector, attached via ethernet port to a Lenovo ThinkStation P620 which has the Aquantia AQC107 NBase-T/IEEE 802.3bz Ethernet Controller. We control this camera via a MATLAB script.

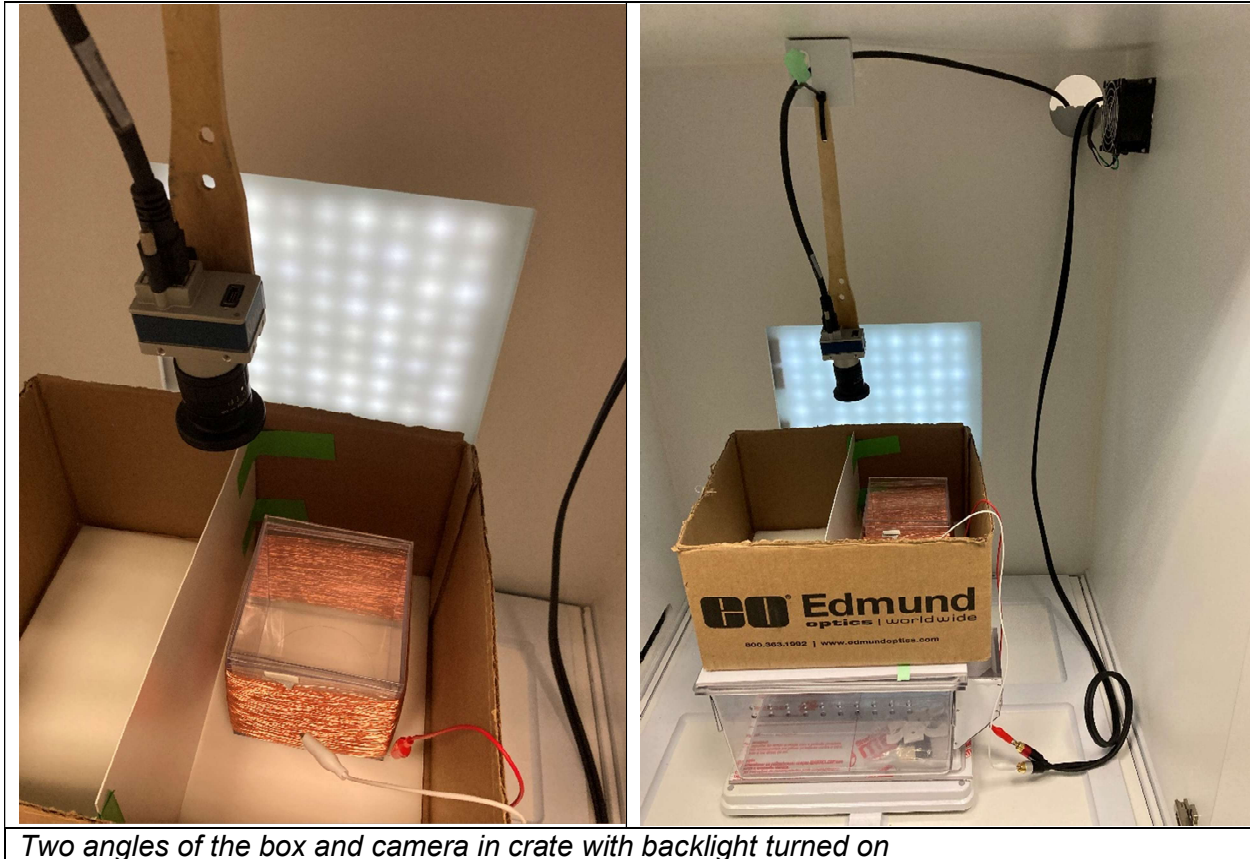

#### Parts List

- Dalsa Genie Nano-M1280 GigE PoE Camera with camera cable
  - <https://www.edmundoptics.com/p/m1280-12-monochrome-dalsa-genie-nano-poe-camera/4048/>
- Lenovo ThinkStation P620
- PoE injector, with ethernet cable
  - <https://www.amazon.com/Single-Port-Injector-Power-Ethernet/dp/B07KFKYF79>
- Crate for video recording (We used the Open Field Test Package w/NIR Lighting for Rat)
  - <https://www.med-associates.com/product/rat-open-field-test-package-with-nir-lighting/>
- Sapera CamExpert image acquisition software

## Assembly

1. Connect the camera to the PoE injector and the PoE to the computer. A popup should appear on the monitor indicating that the camera has connected.
2. Use the MATLAB script described in Section 4 to preview the camera and make recordings.

## Notes

- Make sure to follow the steps in the GigE Vision Quick Start Configuration Guide ([https://www.mathworks.com/matlabcentral/answers/uploaded\\_files/41167/GigEVisionQuickStart.pdf](https://www.mathworks.com/matlabcentral/answers/uploaded_files/41167/GigEVisionQuickStart.pdf)). In specific, ensure that the correct driver from the manufacturer is installed. The driver we used is the Marvell AQtion Windows 64-bit driver (<https://www.marvell.com/support/downloads.html>). In addition, ensure that the firewall is disabled as the firewall may prevent the camera from connecting to the computer. To do so, type in the command line (with administrator controls) `netsh advfirewall set allprofiles state off`.
- Ensure that the light is on when recording.

## Section 4: Computer Control

### Control of the DAQ and Camera via MATLAB

Both the DAQ and camera are controlled by MATLAB scripts.

#### Parts List

- The following MATLAB add-ons are required for control of the DAQ:
  - Signal Processing Toolbox
  - Data Acquisition Toolbox
  - Data Acquisition Toolbox Support Package for National Instruments NI-DAQmx Devices
- The following MATLAB add-ons are required for control of the GigE camera:
  - Image Acquisition Toolbox
  - Image Processing Toolbox
  - Image Acquisition Toolbox Support Package for GigE Vision Hardware
  - Image Acquisition Toolbox Support Package for Teledyne DALSA Sapera Hardware

## Assembly

1. Download and install all required add-ons.
2. To control the DAQ:
  - a. Open `/code/DAQ_code/DAQ_HTR_Code`
  - b. The first section allows the user to enter the number of seconds the session will run for. All succeeding sections (which are adapted from [Revenaga et. al, 2019](#)) apply filters and thresholds to the raw data, correct the raw voltage signal to baseline, and transform the signal to absolute value. As is currently written, the top threshold value is given as 0.2 V, the minimum distance between events is 4000 ms, and the maximum width of an event is 1000 ms. These values worked well for us with few false positives and negatives but could be altered if needed.
  - c. Pasted on the graph will be the number of peaks (head twitches) recorded.
  - d. Pasted in the command window will be the time points of the head twitches.

- e. Save and export the data.mat file. Do not run the script again unless the data.mat file has been saved and exported; otherwise, the data file will be overwritten.
3. To control the camera:
  - a. Within the /code/HTR-recording-master folder, there is a file titled “master\_head\_twitch” that controls both the preview functionality and the recording functionality of the camera. When running this file, make sure that all files within HTR-recording-master folder and the subfolders within are also added to the current MATLAB path. For instance the “bin” folder has two important functions (“newid” and “stoploop”) that are necessary for the program to run.
4. To control the DAQ and camera simultaneously:
  - a. While it is possible to institute a trigger so that the DAQ and camera can run at the exact same time, we have found that is sufficient to open two windows of MATLAB and run the DAQ script soon after running the camera recording script.


#### Notes

- Ensure that all necessary folders and subfolders (such as the “bin” file and “Boubat-UpdateableText-90a229b” file) are added to the correct path.

## Addendum: Scaling Up

### Changes to the Setup when Scaling Up

- It is useful to run multiple mice simultaneously, which can be accomplished by scaling up the setup as shown in the picture below:

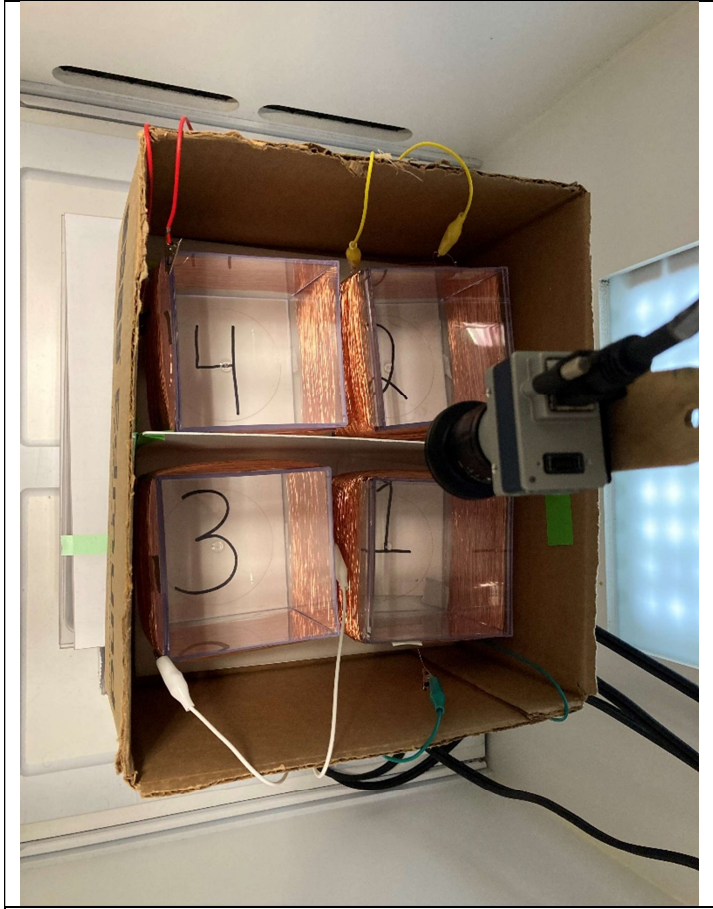

*Setup in which the head twitches of four mice can be automatically quantified simultaneously*

- Similar steps as described in this manual can be applied in order to add more boxes.
- To run the code for two boxes, run `/code/DAQ_code/DAQ_2_Boxes`. Be sure to specify which analog inputs are being used.
- Similarly, run `/code/DAQ_code/DAQ_3_Boxes` or `/code/DAQ_code/DAQ_4_Boxes` to run the code for three or four boxes, respectively.